TENS/EMS Muscle Stimulator Quick Guide

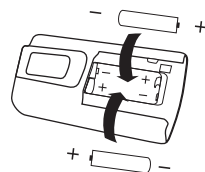
Art. no 34-1548 Modell LT3011A

Please read the entire instruction manual before use and save it for future reference. We reserve the right for any errors in text or images and for making any necessary technical changes to this document. If you should have any questions concerning technical problems please contact our Customer Services.


NOTE! Read all the safety instructions in the User Guide before using the machine. This is only a Quick Guide!

Use


1. Insert the batteries (4 x AAA/LR03) as shown, in the battery compartment.

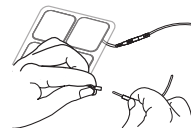2. Connect the cables to the electrode pads.

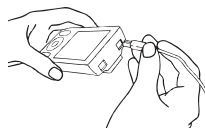
3. Connect the electrode pads' cables to the TENS/EMS machine (first make sure the device is switched off).

4. First wash the skin thoroughly with soap and water, wipe dry and then attach the self-adhesive electrode pads. Make sure they make full contact over the entire surface. **Note!** See section 11. *Placement of the electrode pads* for information in the User Guide for exact placement of the electrode pads on the body.
5. Switch on the machine: Press [⏻], the display lights up. **Warning!** Before switching on the machine, read the warning messages in the User Guide. Pay particular attention to circumstances when the machine must not be used.

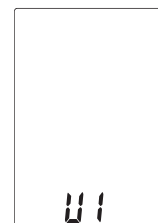6. Select the user by pressing [U1/U2] to select U1 or U2. Press [S] to confirm the selected user. There is a memory for 2 different users, so that 2 people can use the product with their individual settings.

Choose a treatment and area

There are 3 different treatments: **TENS**, **MASSAGE** or **EMS**. Press [M] to select the required treatment.

There are 9 different treatment areas.

Press [S] to select the treatment area:

Shoulder, Neck, Back, Elbow, Hip, Ankle, Foot, Wrist and Knee.

The selected treatment and treatment area flash.

Select the treatment programme

Each treatment area has 3 or 4 programmes for each treatment. The treatment programmes have a different wavelength, pulse rate and pulse width. Diagrams with a description of all programmes can be found at the end of the User Guide.

Press [P] and select the treatment programme (**P1**, **P2**, **P3** or **P4**) for the treatment area you have chosen. The symbol for the selected programme is shown at the bottom of the display.

Start the treatment and adjust the intensity

1. Press [CH1 +] or [CH 2 +] to start the treatment. Press repeatedly to increase the output intensity, each press increases the intensity in 1 mA increments.
2. Press [CH1 -] or [CH2 -] repeatedly to reduce the output intensity for channel 1 or channel 2. Adjustment occurs in 1 mA increments.
3. The set intensity (mA) and the remaining treatment time are shown at the bottom of the display.

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Pausing the treatment

Press [⏻] to pause the treatment. Press again to continue.

Switch off

1. Press and hold [⏻] for about 3 seconds to switch off. **Note!** Except for working mode, the device will automatically turn off when no button is pressed in 3 minutes.
2. Carefully remove the electrode pads from the body and fit the protective backing plastic.
3. Wash the skin with mild soap and water when the electrode pads have been removed.
 - Store the electrode pads in a plastic bag.
 - Seek medical attention if the skin is irritated.

Keypad lock

Note! If you do not press any buttons for 20 seconds this will enable a keypad lock that locks the buttons so that the intensity cannot be increased and no button functions can be used without disabling the keypad lock.

A key symbol is shown on the display when the security lock is enabled.

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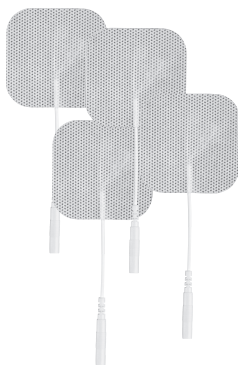
TENS/EMS Muscle Stimulator

Art.no 34-1548 Model LT3011A

Please read the entire User Guide before use and save it for future reference. Errors and omissions excepted. Technical specifications subject to change without notice. If you should have any questions concerning technical issues please contact our Customer Services.

Table of Contents

1. Safety	5
2. The package contains	7
3. Product description	7
4. Display	9
5. Functions	10
<i>5.1 Treatments</i>	<i>10</i>
<i>5.2 Treatment areas</i>	<i>10</i>
<i>5.3 Treatment programmes</i>	<i>10</i>
6. Getting started	10
<i>6.1 Insert the batteries</i>	<i>10</i>
<i>6.2 Connect the cables to the electrode pads</i>	<i>11</i>
<i>6.3 Connect the cables to the TENS/EMS machine</i>	<i>11</i>
<i>6.4 Place the electrode pads on the skin</i>	<i>11</i>
<i>6.5 Switch on the TENS/EMS machine</i>	<i>12</i>
<i>6.6 Select the user</i>	<i>12</i>
<i>6.7 Choose a treatment and treatment area</i>	<i>12</i>
<i>6.8 Select the treatment programme</i>	<i>12</i>
<i>6.9 Start the treatment and adjust the intensity</i>	<i>13</i>
<i>6.10 Pause the treatment</i>	<i>13</i>

6.11 <i>Switch off the machine</i>	13
6.12 <i>Quick stop</i>	13
6.13 <i>Keypad lock</i>	13
6.14 <i>Check what is stored in memory</i>	14
6.15 <i>Battery warning</i>	14
6.16 <i>Individual settings U1</i>	14
7. Care and Maintenance	15
8. Troubleshooting	15
9. Recycling	16
10. Specifications	16
11. Placement of the electrode pads	17
Appendix 1 TENS	75
Appendix 2 MASSAGE	76
Appendix 3 EMS	78
Appendix 4 EMC	79

1. Safety

Warning! None of the product parts may be dismantled or modified in any way. Check all the product parts before each use and replace if a part is damaged in any way. Only original parts may be used!

- The product is intended for home use.
- Follow the user advice in the User Guide. Keep in mind the restrictions and risks associated with this product. Read the warning texts and symbols on the product.
- The product must not be used in the following circumstances:
 - If a heart condition has been established or is suspected.
 - If skin irritation is experienced or hypersensitivity with electrical stimulation or the use of gel.
 - If a rash or persistent pain occurs, stop the treatment and seek medical attention.
- Placement of the electrode pads and settings must be based on advice from a doctor or a person experienced in administering treatment involving similar products.
- Local irritation of the skin may occur with prolonged use of the electrode pads in the same place.
- Do not use the electrode pads on swollen, infected, inflamed or broken skin. Make sure the electrode pads have good contact with the skin.
- If you find the stimulation level uncomfortably high, the impulses must be reduced to a more agreeable level. Stop the treatment and consult a doctor if you are unsure.
- The product may not be used by persons regularly out in traffic or persons who monitor machines and processes.
- Do not use this product in premises where aerosols or oxygen are stored or used.
- Do not use this product in the vicinity of explosive or highly flammable substances or gases.
- Do not use this product at the same time as other products that send electrical impulses to the body.
- Do not connect the electrode pads to any other device or do not connect any other device, i.e. earphones or the like to the TENS/EMS machine.
- Do not use sharp objects to press the buttons.
- Ensure that electrode pads and cables are undamaged before use.
- The TENS/EMS machine may only be used in conjunction with the accompanying accessories.
- Treatment may not be given above the neck, i.e. to the face, on the throat, or anywhere on the head.

Equipment capable of delivering output values in excess of 10 mA r.m.s. or 10 V r.m.s. averaged over any period of 5s.

Applied part of type BF.

Disposal in accordance with Directive 2012/19/EU (WEEE).

Complies with the European Medical Device Directive (93/42/EEC) and amended by directive 2007/47/EC requirements. Notified body TÜV Rheinland (CE0197).

The name and the address of the manufacturer.

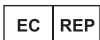
Manufacturer:

Name: Shenzhen Dongdixin Technology Co., Ltd.

Add: No. 3 Building XiliBaimang Xusheng industrial

Estate 518108, Nanshan Shenzhen, China

Tel: 0086-755-27652316 Fax: 0086-755-27652674

The name and the address of the authorized EC-representative.

Authorized EC-representative:

Shanghai International Holding Corp. GmbH (Europe)

Eiffestraße 80, 20537 Hamburg Germany

Tel: 0049-40-2513175 Fax: 0049-40-255726

Batch code

Serial number

Refer to Instruction Manual.

This product must not be used in the following circumstances

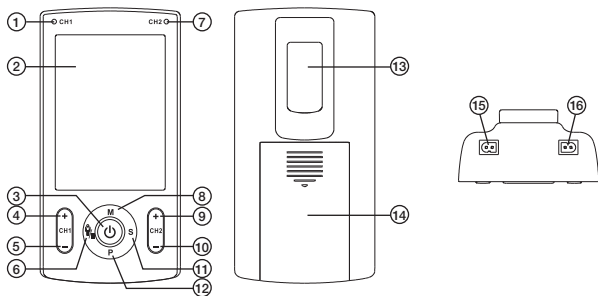
- Pregnant women should consult a doctor before use.
- The product must not be used if the treatment area has been affected by carcinogenic changes.
- The product must not be used on persons who have a pacemaker or defibrillator.
- The product must not be used on persons who have diabetes, abnormally high blood pressure, or who suffer from epilepsy.
- The product must not be used on persons who have serious circulation problems in lower parts of the body, or hernias in the groin or abdominal area.
- The product must not be used if myocarditis has been diagnosed.
- Do not use the product for symptoms that have not been diagnosed.

Warning! Do not use the product if you are unsure how it works. In the event of uncertainty, consult a doctor or a person experienced in administering treatment involving similar products.

2. The package contains


- TENS/EMS machine
- 2 electrode cables, each with 2 connections for electrode pads
- 4 electrode pads, 50 × 50 mm
- User Guide and Quick Guide
- 4 × AAA/LR03 batteries

3. Product description

1. CH1 – LED, lights when channel 1 is enabled
2. LCD display
3.  Switch
 - Press to switch on, the display lights up.
 - Press and hold for about 3 seconds to switch off, the display goes out.
 - In treatment mode: Stop the treatment.
 - In memory mode: Exit memory mode.
 - In setting mode: Exit from the setting state.
4. [CH 1+] Increase the intensity on channel 1. Adjustment from 0 (off) to max.
5. [CH 1-] Decrease the intensity on channel 1. Adjustment from max to 0 (off).
In setting mode: Select the next parameter.
6.  Select user, check data in memory or pause the treatment.
 - Select user U1 (user 1) or U2 (user 2).
 - In memory mode: Check data from previous treatments stored in memory.
 - In treatment mode: Pause and resume the treatment.
7. CH2 – LED, lights when channel 2 is enabled
8. [M] In standby mode, select the treatment: TENS – MASS – EMS.
9. [CH 2+] Increase the intensity on channel 2. Adjustment from 0 (off) to max.
In setting mode: Adjust the user program parameter.
10. [CH 2-] Decrease the intensity on channel 2. Adjustment from max to 0 (off).
In setting mode: Adjust the user program parameter.
11. [S] Confirm the selected users, select the treatment area, open or delete data in memory.
 - Confirm the selected user U1 or U2.
 - Choose the treatment area. Press to step between 9 different treatment areas.
 - Press and hold for 3 seconds to open memory mode (in standby mode).
 - Press and hold for 3 seconds and confirm deletion of memory data (in memory mode).
12. [P] Standby mode: Select the treatment programme P1, P2, P3 or P4 before the treatment starts.
 - The programmes have a different wavelength, pulse rate and pulse width. Diagrams with a description of all programmes can be found at the end of the User Guide.
13. Belt clip
14. Battery cover
15. Connection for channel 1
16. Connection for channel 2

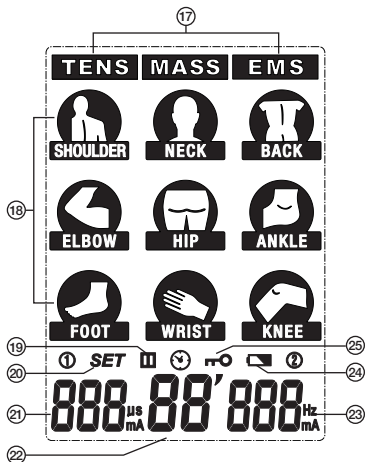4. Display

17. TENS, MASS or EMS shows the chosen treatment.
18. Treatment area on the body: Shoulder, neck, back, elbow, hip, ankle, foot, wrist and knee. (Does NOT show where the electrode pads should be placed on the body).
19. Pause indicator
20. Settings mode symbol
21. Pulse duration and intensity for channel 1.
22. Treatment duration
23. Pulse frequency and intensity for channel 2.
24. Battery charge status indicator.
25. The intensity level keypad lock is activated.

Note! If the electrode pads are not attached to the skin, the intensity can only be adjusted up to 9, at 10 the machine returns to 0 again. When the electrode pads are attached to the skin, the intensity can be adjusted between 0 and maximum value.

5. Functions

5.1 Treatments

Electrical muscle stimulation is used by more and more sportsmen and sportswomen to improve their performance. This is achieved by increasing muscle strength, volume and endurance.

- TENS (Transcutaneous Electrical Nerve Stimulation): Electro therapy with low-intensity and short impulses to e.g. relieve pain.
- EMS (Electrical Muscle Stimulation): Electrical muscle stimulation which increases blood circulation.
- MASSAGE Electrical stimulation that stimulates muscles in the treatment area to mild activity which in turn results in an improved sense of well-being.

5.2 Treatment areas

Shoulder, neck, back, elbow, hip, ankle, foot, wrist and knee.

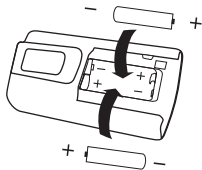
5.3 Treatment programmes

- TENS has 3 programmes: P1, P2 and U1 for each treatment area.
 - EMS has 3 programmes: P1, P2 and U1 for each treatment area.
 - MASSAGE has 4 programmes: P1, P2, P3 and P4 for each treatment area.
- The programmes have a different wavelength, pulse rate and pulse width.

6. Getting started

6.1 Insert the batteries


1. Open the battery cover by pressing it in while sliding in the direction of the arrow.
2. Insert 4 × AAA/LR03 batteries in the battery compartment as shown by the marking.

6.2 Connect the cables to the electrode pads


1. Insert the cable connectors into the electrode pad connectors.
2. Push the connectors together so that no metal is visible on any connection.

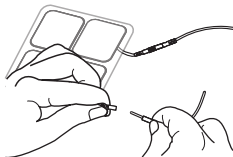Warning! Only use electrode pads that are approved according to: IEC/EN60601-1 and IEC/EN60601-1-2.

6.3 Connect the cables to the TENS/EMS machine

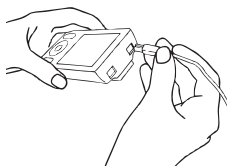1. First make sure the machine is switched off.
2. Grasp an insulated part of the connector and connect this to one of the outlets (15) or (16) on the top.
3. It is possible to use one channel with two electrode pads, or both channels and thus with four electrode pads at the same time. If you use both channels, two different areas of the body can be stimulated simultaneously.

6.4 Place the electrode pads on the skin

The electrode pads should only be used by one person for reasons of hygiene. Additional electrode pads are available as an accessory (our art. no 34-4103-1). Make sure that the skin is completely dry and clean.

1. Remove the protective backing plastic from the self-adhesive surface. Save the protective backing plastic and refit to the electrode pads after use.
2. Wash the skin thoroughly with soap and water, wipe dry and then attach the self-adhesive electrode pads. Make sure they make full contact over the entire surface.
3. Start the TENS/EMS machine when the electrode pads have been attached to the skin.

Warning! Never remove the electrode pads from the skin when the TENS/EMS machine is switched on.

Note! See section 11. *Placement of the electrode pads* for information about the exact placement of the electrode pads on the body.

6.5 Switch on the TENS/EMS machine

Read the safety instructions before switching on the machine. Pay particular attention to circumstances when the machine must not be used. The TENS/EMS machine is not a toy, but a powerful treatment aid.

Switch on the machine: Press , the display lights up.

6.6 Select the user

There is a memory for 2 different users (U1) and (U2), so that 2 people can use the product with individual settings.

1. Press  to select user 1 or 2.
2. Press [S] to confirm the selected user.

6.7 Choose a treatment and treatment area

1. There are 3 different treatments: **TENS**, **MASSAGE** or **EMS**. Press [M] to select the required treatment.
2. There are 9 different treatment areas. Press [S] to select the treatment area.
3. The selected treatment and the treatment area flash.

6.8 Select the treatment programme

For the 3 treatments and the 9 treatment areas there is a total of 78 treatment programmes. TENS and EMS has 27 programmes and MASSAGE has 24 programmes.

Press [P] and select the treatment programme (**P1**, **P2**, **P3**, **U1*** or **P4**) for the treatment area you have chosen. *P1*, *P2*, *P3*, *U1* or *P4* is shown at the bottom of the display for the selected programme.

The treatment programmes have a different wavelength, pulse rate and pulse width. Diagrams with a description of all programmes can be found at the end of the User Guide.

* See section 6.16 *Individual settings U1*

6.9 Start the treatment and adjust the intensity

1. Press [CH1 +] or [CH 2 +] to increase the output intensity from **0** and upward for channel 1 or channel 2 and in doing so, start the treatment.
2. Press [CH1 -] or [CH2 -] repeatedly to reduce the output intensity for channel 1 or channel 2. Adjustment occurs in 1 mA increments.
3. The set intensity (mA) and the remaining treatment time are shown at the bottom of the display.

The maximum treatment duration is 90 minutes.

Warning! Reduce the stimulation level if you experience discomfort. Seek medical attention if the discomfort persists.

6.10 Pause the treatment

Press  to pause the treatment. Press again to continue.

6.11 Switch off the machine

1. Press and hold  for about 3 seconds to switch off. **Note!** Except for working mode, the device will automatically turn off when no button is pressed in 3 minutes.
2. Carefully remove the electrode pads from the body and fit the protective backing plastic. If the treatment is to be continued within a short time, you may leave the electrode pads on the body.
3. Wash the skin with mild soap and water when the electrode pads have been removed.

Store the electrode pads in a plastic bag.

Seek medical attention if the skin is irritated.

6.12 Quick stop

If you need to switch off the machine quickly:

Press  to pause the treatment, press again if you wish to continue.

Note!

- The strength of stimulation must be adjusted individually for each user.
- The intensity can be adjusted from **0** (off) to max.
- Switch off the machine before the electrode pads are removed from the body.
- The intensity goes down to 0 automatically when the treatment is finished.

6.13 Keypad lock

Note! If you do not press any buttons for 20 seconds this will enable a security lock/keypad lock that locks the buttons so that the intensity cannot be increased and no button functions can be used without disabling the keypad lock. A key symbol is shown on the display when the keypad lock is enabled. To turn off the keypad lock: Press [CH1-] or [CH2-].

6.14 Check what is stored in memory

The machine saves different parameters in memory such as the memory number of the treatment and its programme, treated body part, average intensity and duration. A maximum of 30 treatments may be stored in memory, after 30 treatments, the oldest is deleted and the most recent saved.

To check what is stored in memory: Press and hold [S] for 3 seconds to open memory mode. Press [🔍/■] to read the various memory parameters. Press [CH1] to go to the next memory number.

Note!

- The machine returns to standby mode when you press [🔌], or automatically after 30 seconds.
- Press and hold [S] for 3 seconds to show the total number of treatments in the memory, they are shown flashing on the display:
 - Delete everything in memory: press [🔍/■].
 - Cancel deletion and return to standby mode: press [🔌].

6.15 Battery warning

When the battery symbol flashes, the batteries should be replaced as soon as possible. The machine will continue to work for a few hours before shutting down.

6.16 Individual settings U1

In TENS and EMS modes there is a mode (U1) for programming your own individual treatment settings.

1. Turn the machine on and select U1 by pressing [P].
2. Keep [P] held in to enter the settings mode for U1.
3. Select the desired parameters using [CH1-].
4. Adjust the selected parameter using [CH2+] or [CH2-].
 - Adjust the pulse duration from 50 to 350 μ s.
 - Adjust the pulse frequency from 1 to 150 Hz.
 - Adjust the treatment duration from 5 to 90 minutes.
5. Press [🔌] to confirm settings and enter standby mode.

Note: Reset U1 to default factory settings by holding in [CH2-] at the same time as you insert the batteries.

7. Care and Maintenance

Remove the batteries before cleaning. Wipe off the machine with a soft, damp cloth. Use a mild detergent if necessary, never a strong detergent or solvent.

Note!

- Check the electrode pads before each use, and replace if damaged in any way.
- Store the electrode pads in a plastic bag between treatment sessions.

8. Troubleshooting

Problem	Probable cause	Solution
The display does not light up.	The batteries are spent or inserted incorrectly.	<ol style="list-style-type: none"> 1. Replace the batteries. 2. Check whether the batteries have been installed correctly.
Weak stimulation.	<ol style="list-style-type: none"> 1. The electrode pads have dried or are dirty. 2. Wrong placement. 3. Faulty cables. 	<ol style="list-style-type: none"> 1. Replace with new electrode pads. 2. Check the electrode pad placement. 3. Check the cables and ensure they are connected properly.
Stimulation is uncomfortable.	<ol style="list-style-type: none"> 1. Intensity too high. 2. The electrode pads sit too close together. 3. Damaged electrical pads or cables. 4. The electrode pads' active area is too small. 	<ol style="list-style-type: none"> 1. Reduce the intensity. 2. Change the placement. 3. Replace them. 4. Replace the electrode pads if smaller than 4 × 4 cm.
Intermittent output.	<ol style="list-style-type: none"> 1. Faulty cables or connection. 2. Program function. 	<ol style="list-style-type: none"> 1. Check the cable connections, push them in fully. Reduce the intensity. Disconnect the cable and reconnect. Replace cables that are faulty. 2. Some programmes may be perceived as intermittent, read the programme description.
Stimulation is ineffective.	Placement of the electrode pads incorrect, check.	Read the description and change the placement of the electrode pads.

9. Responsible disposal

This symbol indicates that this product should not be disposed of with general household waste. This applies throughout the entire EU. In order to prevent any harm to the environment or health hazards caused by incorrect waste disposal, the product must be handed in for recycling so that the material can be disposed of in a responsible manner. When recycling your product, take it to your local collection facility or contact the place of purchase. They will ensure that the product is disposed of in an environmentally sound manner.

10. Specifications

Material	ABS plastic
Channels	2 channels
Power supply	4 × AAA/LR03 batteries
Waveform	Two phases square wave
Pulse length	30–350 µs
Pulse frequency	1–290 Hz
Treatment duration	5–90 minutes
Intensity	Adjustable from 0 to 90 mA (at 1000 ohms)
Operating environment	5 to 40 °C at 30–85% RH and air pressure 700–1060 hPa
Storage	-10 to 50 °C at 10–90% RH and air pressure 700–1060 hPa
Dimensions	117 × 60 × 25 mm (without belt clip)
Weight	110 g (without batteries)

Warning! Read the entire User Guide before use.

Unit type BF.






The product must be handled in accordance with decree 2012/19/EU – WEEE (Waste Electrical and Electronic Equipment).





CE 0197

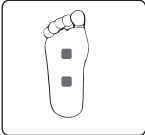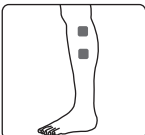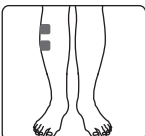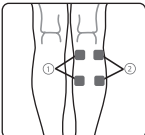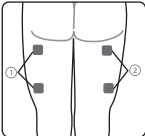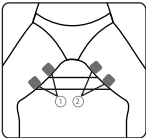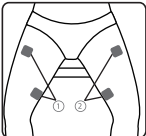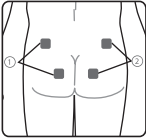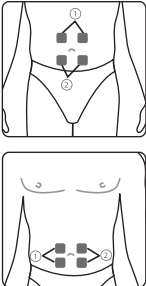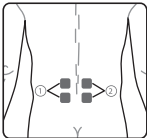
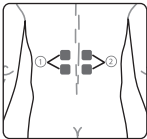
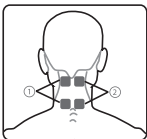
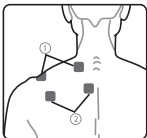
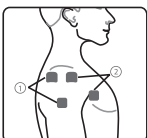
The product is in conformity with the guidelines set down in the European Medical Device Directive (93/42/EEC) and as amended by Directive 2007/47/EC. Notified Body TÜV Rheinland (CE0197).

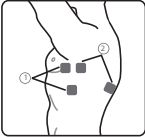
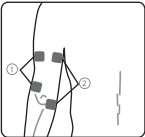
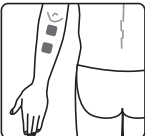
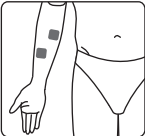
The product conforms to the requirements set out in the following European standards: EN 60601-1, EN 60601-1-2, EN 60601-2-10, ISO10993-5, ISO10993-10

11. Placement of the electrode pads

Electrode placement	Treatment position
	<p>In upright sitting position with feet resting on the floor.</p>
	<p>Upright sitting position with feet resting on the floor.</p>
	<p>Upright sitting position with feet under a piece of furniture to avoid bending of ankles.</p>
	<p>Upright sitting position with support for feet and back (e.g., sitting in a door frame).</p>
	<p>Lying flat on one's stomach with ankles comfortably immobilised.</p>

	<p>Upright sitting position, knees pressed together around a hard (but comfortable) object.</p>
	<p>Upright sitting position, knees pressed together around a hard (but comfortable) object.</p>
	<p>Lying on one's stomach or standing up.</p>
	<p>Lying on one's stomach or standing up.</p>

	Upright sitting position.
	Upright sitting position.
	Upright sitting position.
	Upright sitting position.
	Upright sitting position with elbows pressed against side of the body (inside chair armrests).

	<p>Upright sitting position with elbows pressed away from the body (outside chair armrests).</p>
	<p>Upright sitting position with arms resting on chair armrests.</p>
	<p>Upright sitting position with lower arms resting on chair armrest and palms of hands upward.</p>
	<p>Upright sitting position with lower arms resting on chair armrests and palms resting on chair armrests.</p>
	<p>Upright sitting position with lower arms resting on chair armrests and palms of hands upward. Hands clasping some object so fingers are slightly bent.</p>

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Appendix/Bilaga/Liite 1 TENS

Program Therapy part	P1			P2			U1 (Default)			Treatment time (min)
	Pulse rate (Hz)	Pulse width (µs)	Waveform	Pulse rate (Hz)	Pulse width (µs)	Waveform	Pulse rate (Hz)	Pulse width (µs)	Waveform	
SHOULDER	2~80	100~200	Simple modulated pulse	50	200	Synchronous	2~125	100~200	Simple modulated pulse	30
NECK	2	180	Continuous	80	70~180	Pulse width modulated	80	180	Continuous	
BACK	80/2	180	Hans	80	70/180	Hans	100	330/200	Amplitude modulated	
ELBOW	2	180	Continuous	80	70~180	Pulse width modulated	2~125	100~200	Simple modulated pulse	
HIP	100	330/200	Amplitude modulated	125	330/200	Amplitude modulated	80	330/200	Amplitude modulated	
ANKLE	2~60	100~200	Simple modulated pulse	2~8	300	Pulse rate modulated	2~40	100~200	Simple modulated pulse	
FOOT	80	70~180	Pulse width modulated	80	70/180	Hans	2~100	100~200	Simple modulated pulse	
WRIST	50	200	Synchronous	65	200	Synchronous	2~80	100~200	Simple modulated pulse	
KNEE	50	350	Asynchronous	2	180	Continuous	80	200	Alternate Ramped Burst	

Appendix/Bilaga/Liite 2 MASSAGE

Therapy part	program
SHOULDER	P1, P2, P3, P4
NECK	P1, P2, P3
BACK	P1, P2, P3
ELBOW	P1, P2
HIP	P1, P2, P3, P4
ANKLE	P1, P2
FOOT	P1, P2
WRIST	P2, P3
KNEE	P1, P2

Program	Phase	Pulse rate (Hz)	Pulse width (μ s)	Work time (s)	Rest time (s)	Treatment Time(min)
P1	1	8	100	4	/	30
P2	1	25	30~220	3.5	1.0	
	2	25	30~220	2.5	0.9	
	3	33	30~220	1.9	0.9	
	4	43	30~220	1.3	0.8	
	5	53	200	0.9	0.7	
	6	69	200	0.7	0.6	
	7	79	200	0.5	0.5	
	8	69	200	0.7	0.6	
	9	53	200	0.7	0.6	
	10	43	30~220	1.3	0.8	
	11	33	30~220	1.9	0.9	
	12	25	30~220	2.5	0.9	
	13	25	30~220	3.5	1	
P3	1	83	50~220	4	1	
	2	100	50~220	3.8	0.8	
	3	111	50~220	3.1	0.72	
	4	118	50~220	2.6	0.6	
	5	132	50~220	2.3	0.6	

Program	Phase	Pulse rate (Hz)	Pulse width (μ s)	Work time (s)	Rest time (s)	Treatment Time(min)
P3	6	118	50~220	2.6	0.6	30
	7	111	50~220	2.8	0.7	
	8	100	50~220	3.3	0.8	
P4	1	147	30~220~150	12.0	1.0	
	2	169	30~220~150	10.3	0.9	
	3	196	30~220~150	8.5	0.6	
	4	237	30~220~150	6.8	0.6	
	5	285	30~220~150	5.1	0.4	
	6	290	30~220~150	5.7	0.5	
	7	238	30~220~150	6.3	0.5	
	8	197	30~220~150	8	0.6	
	9	191	30~220~150	8.5	0.7	
	10	168	30~220~150	9.1	0.8	
11	150	30~220~150	10.8	0.9		

Appendix/Bilaga/Liite 3 EMS

Program Therapy part	P1			P2			U1 (Default)			Treatment time (min)
	Pulse rate (Hz)	Pulse width (µs)	Waveform	Pulse rate (Hz)	Pulse width (µs)	Waveform	Pulse rate (Hz)	Pulse width (µs)	Waveform	
SHOULDER	80	150	Synchronous	50	200	Synchronous	40	200	Synchronous	30
NECK	1	50	Continuous	2~60	100~200	Simple modulated pulse	2~8	300	Pulse rate modulated	
BACK	50	350	Synchronous	65	350	Synchronous	80	150	Synchronous	
ELBOW	50	150	Synchronous	50	200	Synchronous	50	350	Asynchronous	
HIP	65	350	Synchronous	40	200	Synchronous	50	350	Synchronous	
ANKLE	50	200	Synchronous	65	200	Synchronous	50	200	Synchronous	
FOOT	50	200	Synchronous	50	200	Synchronous	65	200	Synchronous	
WRIST	50	200	Asynchronous	50/8	200	Asynchronous	50	350	Synchronous	
KNEE	50	350	Synchronous	50	350	Synchronous	40	350	Synchronous	

Appendix/Bilaga/Liite 4 EMC

IMPORTANT INFORMATION REGARDING ELECTROMAGNETIC COMPATIBILITY (EMC)

MEDICAL ELECTRICAL EQUIPMENT needs special precautions regarding EMC. For detailed description of EMC requirements please contact Consumer Relations.

Portable and mobile RF communications equipment can affect MEDICAL ELECTRICAL EQUIPMENT.

Non-Life Support Equipment Separation Distance Calculation (3 Vrms / 3 V/m compliance)

Rated maximum output power of transmitter (W)	Separation distance according to frequency of transmitter (m)		
	150 KHz to 80 MHz	80 MHz to 800 MHz	800 MHz to 2.5 GHz
	$d = \left[\frac{3,5}{V_1} \right] \sqrt{P}$	$d = \left[\frac{3,5}{E_1} \right] \sqrt{P}$	$d = \left[\frac{7}{E_1} \right] \sqrt{P}$
0.01	0.12	0.12	0.23
0.1	0.37	0.37	0.74
1	1.17	1.17	2.33
10	3.69	3.69	7.38
100	11.67	11.67	23.33

Guidance and manufacturer's declaration — electromagnetic immunity

The device is intended for use in the electromagnetic environment specified below. The customer or the user should assure that it is used in such an environment.

Immunity test	IEC 60601 test level	Compliance level	Electromagnetic environment - guidance
Electrostatic discharge (ESD) IEC 61000-4-2	±6 kV contact ±8 kV air	Complies	Floors should be wood, concrete or ceramic tile. If floors are covered with synthetic material, the relative humidity should be at least 30 %.
Radiated RF IEC 61000-4-3	3 V/m 80 MHz to 2.5 GHz	Complies	Field strengths outside the shielded location from fixed RF transmitters, as determined by an electromagnetic site survey, should be less than 3 V/m. Interference may occur in the vicinity of equipment marked with the following symbol: 
Conducted RF IEC 61000-4-6	3 Vrms 150 kHz to 80 MHz	Complies	Separation distance calculation provided below. If a known transmitter is present the specific distance can be calculated using the equations.

Electrical fast transient/burst IEC 61000-4-4	±2 kV for power supply lines ±1 kV for input/output lines	Not applicable	The ME equipment is solely battery powered.
Surge IEC 61000-4-5	±1 kV differential mode ±2 kV common mode	Not applicable	
Power frequency (50/60Hz) magnetic field IEC 61000-4-8	3 A/m	Complies	Power frequency magnetic fields should be at levels characteristic of a typical location in a typical commercial or hospital environment.
Voltage dips, short interruptions and voltage variations on power supply input lines IEC 61000-4-11	>95% dip 0.5 cycle 60% dip 5 cycles 70% dip 25 cycles 95% dip 5 sec.	Not applicable	The ME equipment is solely battery powered.

Guidance and manufacturer's declaration - electromagnetic emissions

The device is intended for use in the electromagnetic environment specified below. The customer or the user should assure that it is used in such an environment.

Emissions test	Compliance	Electromagnetic environment - guidance
RF emissions CISPR 11	Group 1	The device uses RF energy only for its internal function. Therefore, its RF emissions are very low and are not likely to cause any interference in nearby electronic equipment.
RF emissions CISPR11	Class B	Complies
Harmonic emissions IEC 61000-3-2	Not applicable	The ME equipment is solely battery powered.
Voltage fluctuations / flicker emissions IEC 61000-3-3	Not applicable	

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