Electric Toothbrush

Art.no 34-2139-1/-2

Please read the entire instruction manual before using the product and then save it for future reference. We reserve the right for any errors in text or images and any necessary changes made to technical data. If you have any questions regarding technical problems please contact our Customer Services.

Safety


- Only use the provided accessories.
- The toothbrush must not be immersed in water. If this should occur despite all precautionary measures, let the toothbrush stand vertically until the battery compartment drains thoroughly dry before attempting to use it again. Do not attempt to reuse the toothbrush until you are sure it is thoroughly dry.
- The toothbrush is not intended for use by children under 3 years of age. Children using this appliance should be supervised.
- The toothbrush must not be exposed to extreme temperature changes.
- Never use the toothbrush if it has been damaged in any way.
- Never attempt to dismantle the toothbrush. Repairs should only be carried out by a qualified tradesman.
- Never let children play with the battery.

Inserting the battery

1. Remove the battery cover.
2. Insert 1 x AA/LR6 battery (sold separately) with the + pole inserted first.
3. Refit the battery cover.

Note: The battery cover only fits in one direction.

If the toothbrush is not to be used for an extended period, remove the battery in order to avoid any possible risk of leakage.

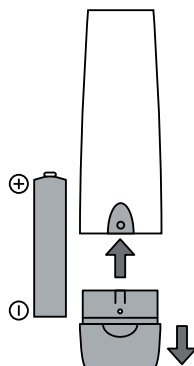
Operating instructions

- Do not use excessive pressure while using the toothbrush. Consult a competent specialist if you are at all uncertain.
- The 2-minute hourglass is meant to be used as a gauge for the optimal amount of time needed for adequate brushing.

Changing the brush head

1. Hold the toothbrush with one hand.
2. Use the other hand to press down on the brush head and twist it anticlockwise to remove it from the handle.
3. Replace with a new brush head by pressing down onto the toothbrush and twisting clockwise.

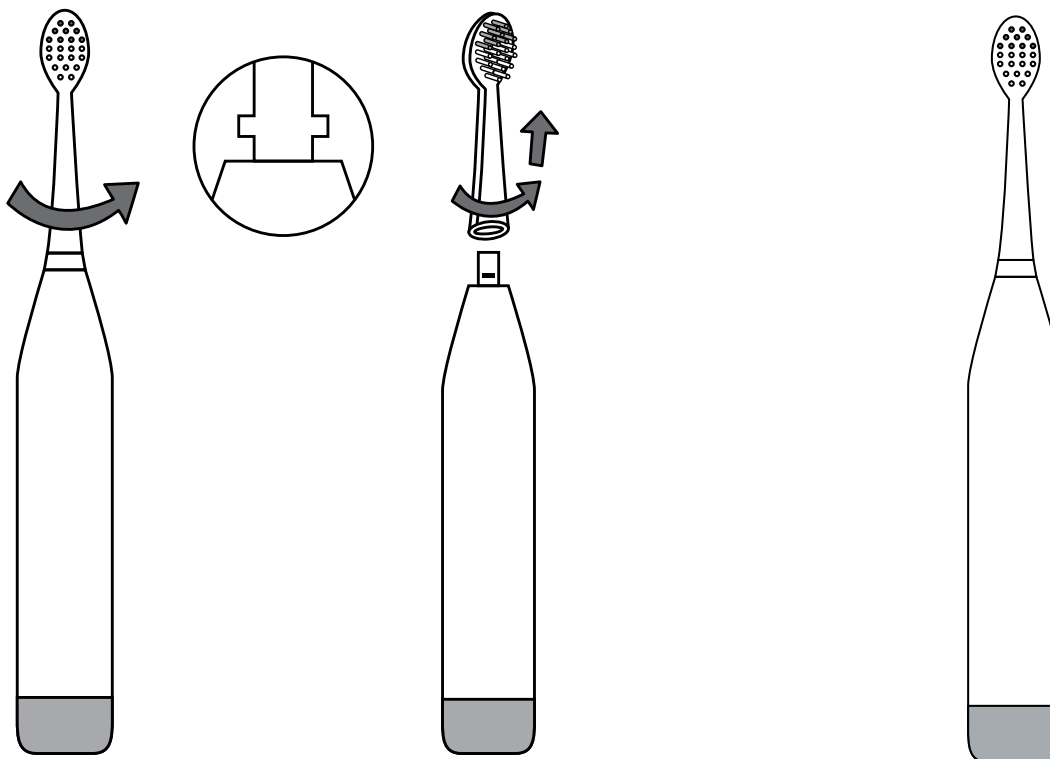Care and maintenance

Rinse the brush head in running water. Remove the battery from the toothbrush if it is not to be used for an extended period.

Disposal

This product should be disposed of in accordance with local regulations. If you are unsure how to proceed, contact your local authority.

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