Travel Toothbrush


Art.no 34-3080-1, -2, -3, -4, -5, -6, -7

Please read the entire instruction manual before using the product and then save it for future reference. We reserve the right for any errors in text or images and any necessary changes made to technical data. If you have any questions regarding technical problems please contact our Customer Services.

Operating instructions


The vibration function is switched on and off by pressing [⏻].

Brush head replacement

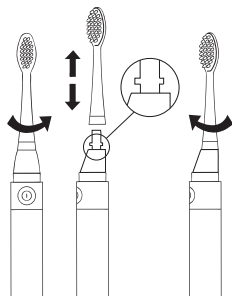1. Twist the brush head anticlockwise and pull it straight upwards to remove it from the handle.
2. Place the new brush head onto the neck of the toothbrush and twist it clockwise until the head is firmly attached.

Battery replacement

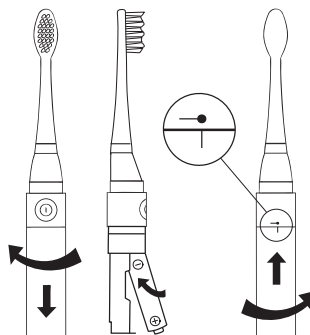1. Twist the handle anticlockwise and remove the housing to get to the battery.
2. Insert a new AAA/LR03 battery.
3. Refit the housing and twist it clockwise to lock it.

Care and maintenance

1. Rinse the brush head in cold water after each use.
2. Shake off the excess water.
3. Take the battery out if the toothbrush is not to be used for a long time.
4. Clean the product by wiping it with a soft, damp cloth.
5. Never immerse or rinse the handle/battery housing in water.

Disposal

This product should be disposed of in accordance with local regulations. If you are unsure how to proceed, contact your local authority.

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