
34-7283

Bathroom Scale

Personvåg • Personvekt • Henkilövaaka

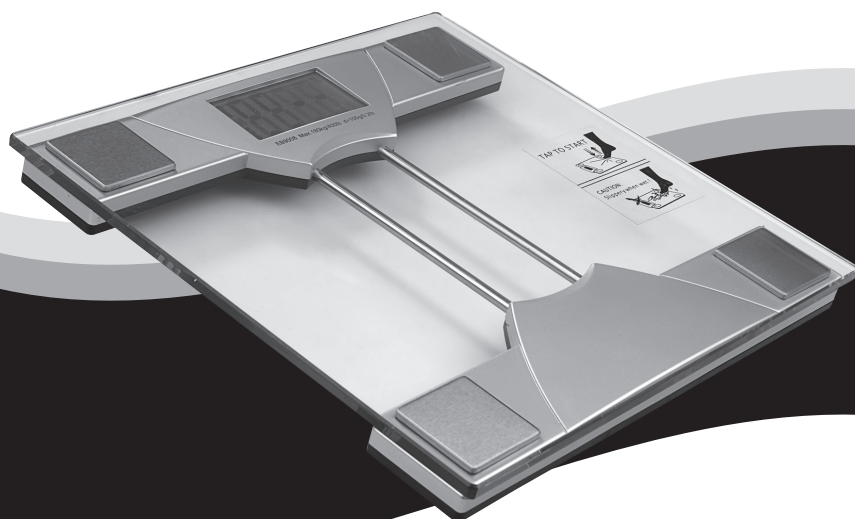Model EB9008

English	3
Svenska	7
Norsk	11
Suomi	15

Bathroom Scale

Article number: 34-7283
Model: EB9008

Please read the entire instruction manual before using and save it for future use.

We apologise for any text or photo errors and any changes of technical data.

If you have any questions concerning technical problems please contact our Customer Service Department (see address on reverse.)

Product Description



- Memory for previous weight records
- Capacity: 180 kg
- Sensitivity: 100 g
- Activate the scale by a light touch of the foot
- Automatic reset to zero
- Automatic shutdown
- Overload indicator

Battery

The battery pack underneath holds 1x CR2032 (3 V) battery (+ pole facing the battery cover). Remove the isolating plastic tag between the battery and the battery holder before first use.


How to change battery: Carefully push the lock tag outward so that the battery loosens, push a new battery (CR2032) down into the holder.

When the battery needs changing "Lo" appears on the display.

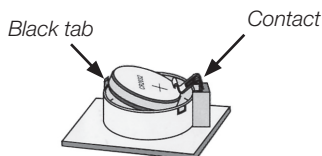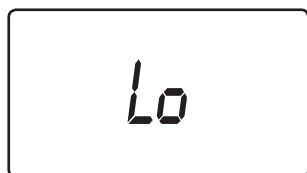Display


If the display shows strange symbols or if the scale does not shut itself off, try to reset it to factory presetting: Remove the batteries and wait at least 3 seconds and then put the batteries back.

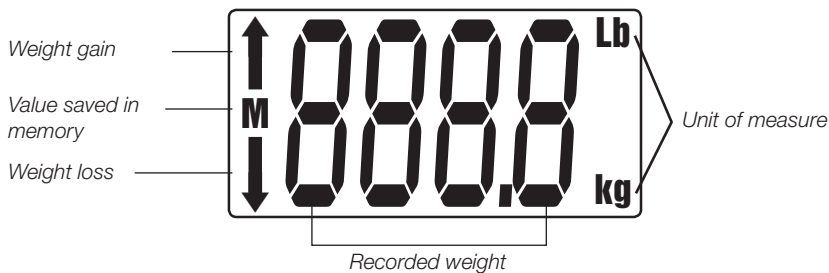Use

Warning! Bare in mind that the glass is slippery when wet, make sure that it is dry before standing on the scale.

1. Place the scale on a firm, even surface (avoid placing the scale on a mat or other soft surface).

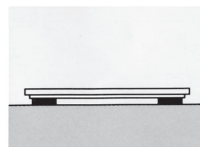
2. Activate the scale by using a little light foot pressure, after a few seconds "0.0" will appear on the display, the scale is now ready for use.


3. Carefully mount the scale, with feet evenly spaced, weight evenly distributed between both legs and between heel and toes. Remain stationary until the display shows a steady reading during a few seconds. The display blinks when reading is done.

4. The read weight is automatically saved in the memory when the display blinks, the memory can save up to 80 weight records.

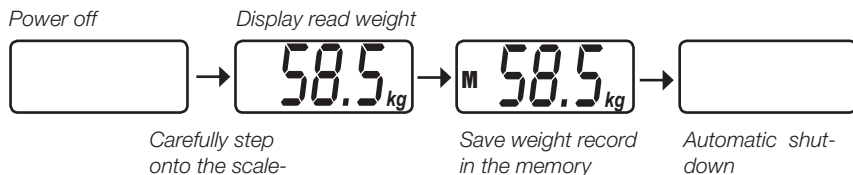
5. The scale have a selector on the bottom with a choice of : Kg-Lb.

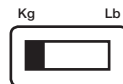Display in KG or Lb

6. The scale compares the most recent record with previous saved records and shows with an arrow up/down if the weight is greater/lesser then the last read weight or the most previous.

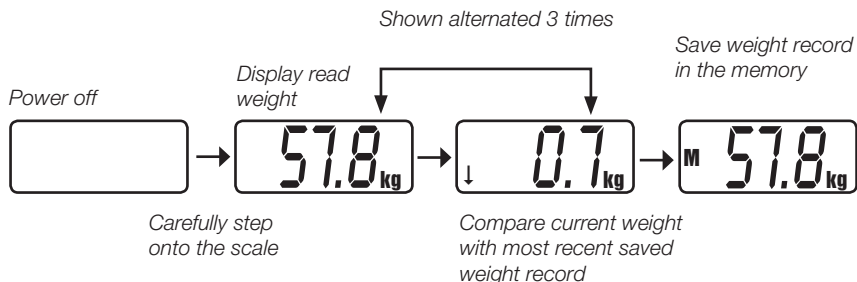
7. If the scale is overloaded "ERR" will appear on the display.

Indication that the weighing failed

If problems occur during the weighing so that the weighing is not shown correctly a "C" will show on the display. Redo the weighing.

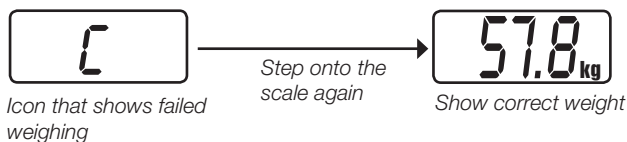Automatic shutdown

The scale automatically shuts off when you step off of it.

Helpful suggestions!

- This product is a precision instrument and should be handled with care and protected from blows and rough treatment.
- Place the scale in the same place every time you weigh yourself.
- The scale need to be placed on a hard, even surface for accurate reading, uneven surfaces help create inaccurate readings.
- Weigh yourself without clothing at the same time of day, preferably in the morning before breakfast. Place your feet in the same place every time you weigh yourself.
- The display rounds up/down to the nearest 100 g.
- Clean the scale using a damp cloth as and when required don't let water run into the scale. Never use strong solvents or detergents.
- Keep the scale away from water and splashes, the electronics can be damaged.
- Do not store the scale on it's side because the display will be on all the time and the batteries will soon discharge.
- Take the batteries out if the scale is not to be used for a long time.
- The scale can be affected by changes in temperature and by being moved. Inaccurate readings may occur if the scale is too cold when used. Allow the scale to stand and acclimatise for approx. 24 hours before using it for the first time.

Protect the environment

Sort your waste and recycle products instead of throwing them away. Electrical and electronic products marked with the recycling symbol should be handled according to local environmental regulations.

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