TENS/EMS Muscle Stimulator Quick Guide


Art.no 34-1548 Modell LT3011A

Please read the entire instruction manual before use and save it for future reference. We reserve the right for any errors in text or images and for making any necessary technical changes to this document. If you should have any questions concerning technical problems please contact our Customer Services.

Note: Read all the safety instructions in the User Guide before using the machine. This is only a Quick Guide!

Use


1. Insert the batteries (4 x AAA/LR03) as shown, in the battery compartment.

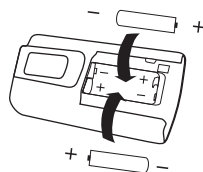2. Connect the cables to the electrode pads.

3. Connect the electrode pads' cables to the TENS/EMS machine (first make sure the device is switched off).

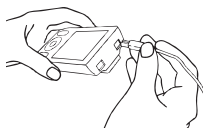
4. First wash the skin thoroughly with soap and water, wipe dry and then attach the self-adhesive electrode pads. Make sure they make full contact over the entire surface.

Note: See section 11. *Placement of the electrode pads* for information in the User Guide for exact placement of the electrode pads on the body.
5. Switch on the machine: Press [⏻], the display lights up.

Warning: Before switching on the machine, read the warning messages in the User Guide. Pay particular attention to circumstances when the machine must not be used.

6. Select the user by pressing [👤/👤] to select U1 or U2. Press [S] to confirm the selected user. There is a memory for 2 different users, so that 2 people can use the product with their individual settings.

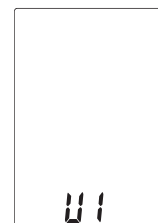Choose a treatment and area

There are 3 different treatments: **TENS**, **MASSAGE** or **EMS**. Press [M] to select the required treatment.

There are 9 different treatment areas. Press [S] to select the treatment area: **Shoulder, Neck, Back, Elbow, Hip, Ankle, Foot, Wrist and Knee.**

The selected treatment and treatment area flash.

Select the treatment programme

Each treatment area has 3 or 4 programmes for each treatment. The treatment programmes have a different wavelength, pulse rate and pulse width. Diagrams with a description of all programmes can be found at the end of the User Guide.

Press [P] and select the treatment programme (**P1**, **P2**, **P3** or **P4**) for the treatment area you have chosen. The symbol for the selected programme is shown at the bottom of the display.

Start the treatment and adjust the intensity

1. Press [CH1 +] or [CH 2 +] to start the treatment. Press repeatedly to increase the output intensity, each press increases the intensity in 1 mA increments.
2. Press [CH1 -] or [CH2 -] repeatedly to reduce the output intensity for channel 1 or channel 2. Adjustment occurs in 1 mA increments.
3. The set intensity (mA) and the remaining treatment time are shown at the bottom of the display.

Warning: Reduce the stimulation level if you experience discomfort. Seek medical attention if the discomfort persists.

Pausing the treatment

Press [👤/👤] to pause the treatment. Press again to continue.

Switch off

1. Press and hold [⏻] for about 3 seconds to switch off.

Note: Except for working mode, the device will automatically turn off when no button is pressed in 3 minutes.
2. Carefully remove the electrode pads from the body and fit the protective backing plastic.
3. Wash the skin with mild soap and water when the electrode pads have been removed.
 - Store the electrode pads in a plastic bag.
 - Seek medical attention if the skin is irritated.

Keypad lock

Note: If you do not press any buttons for 20 seconds this will enable a keypad lock that locks the buttons so that the intensity cannot be increased and no button functions can be used without disabling the keypad lock.

A key symbol is shown on the display when the security lock is enabled.

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