Bathroom Scale

Personvåg • Personvekt • Henkilövaaka

coline®

Art. No.
34-8394

Model
BS-801

| | |
|---------|----|
| English | 3 |
| Svenska | 6 |
| Norsk | 9 |
| Suomi | 12 |

Bathroom Scale

Article number: 34-8394 Model: BS-801

Please read the entire instruction manual before using the product and save it for future reference. We reserve the right for any errors in text or images and any necessary changes made to technical data. If you have any questions regarding technical problems please contact Customer Services (see address on reverse).

1. Product description


- Capacity: 180 kg / 396 lb
- Sensitivity: 100 g / 0.2 lb
- The scale is activated when weight is put upon it.
- Automatic zeroing
- Automatic shutdown

2. Inserting batteries

When the batteries need changing “Lo” appears on the display.

Batteries: 4 x AAA/R03 batteries.

There are two battery compartments on the bottom with room for 2 x AAA/R03 batteries in each battery compartment. Insert the batteries according to the marking in the battery compartment.

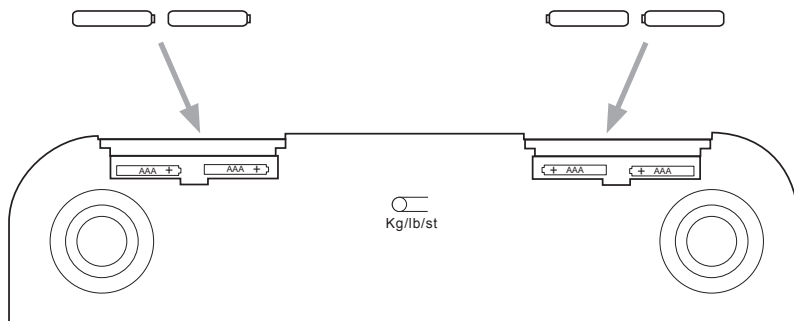3. Display

The display shows the load weight and selected weight units.


Reset

If the display behaves abnormally or if the scale does not shut itself off it may need to be reset to the original factory settings: Remove the batteries and wait at least a few seconds before replacing them.

4. Select weighing units

The scale has a selector button on the bottom that can be used to select weighing units in: kg – st – lb. Each press of the button selects the next weighing unit.

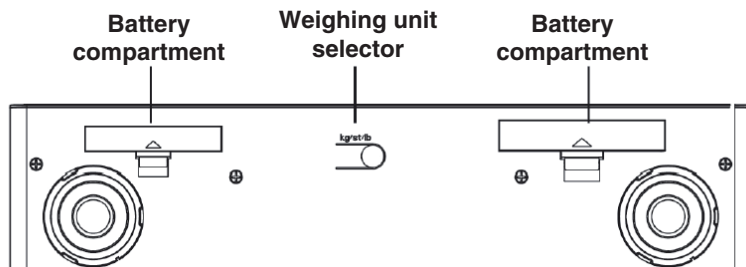5. Operation

1. Place the scale on a firm, even surface (avoid placing the scale on a mat or other soft surface).
2. Carefully mount the scale with your weight evenly distributed between both feet. Remain stationary until the display shows a steady reading for a few seconds. The display will flash once when weighing is complete.

Automatic shutdown

The scale automatically shuts off 10 seconds after you step off it.

7. Troubleshooting

Problem: The display shows ----

Solution: The scale is overloaded. Reduce the weight on the scale.

Problem: "Lo" is shown on the display.

Solution: The batteries are weak. Insert new batteries.

8. Note!

- Keep the scale clear of objects when it is not in use. If the scale has been moved you must allow the display to turn off before mounting the scale. This allows the scale to zero itself.
- This product is a precision instrument and should be handled with care and protected from blows and rough treatment.
- The scale needs to be placed on a hard, even surface for accurate reading. Uneven surfaces help create inaccurate readings.
- Weigh yourself without clothing at the same time of day, preferably in the morning before breakfast. Place your feet in the same place every time you weigh yourself.
- The display rounds up/down to the nearest 100 g.
- Clean the scale using a damp cloth when required. Make sure that no water gets into the scale. Never use strong solvents or detergents.
- Keep the scale away from water and splashes, the electronic components can be damaged.
- Remove the batteries from the scale if it is not to be used for an extended period.
- The scale can be affected by changes in temperature and by being moved. Inaccurate readings may occur if the scale is too cold when used. Allow the scale to stand and acclimatise for approx. 24 hours before using it for the first time.

9. Protecting the environment

Dispose of your waste responsibly, sort and recycle where possible. Electrical and electronic products marked with the recycling symbol should be dealt with according to local environmental regulations.

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BREV Clas Ohlson AB, 793 85 INSJÖN
OSOITE Clas Ohlson Oy, Yrjönkatu 23 A, 00100 HELSINKI